


Clinical Software Innovations thanks you for your interest in **Limb Volumes Professional** (LVP) software. The latest version is LVP 6.0, which offers a number of added capabilities to the already extensive features of LVP5.0. We believe that this software will provide you with a useful and easy to use tool to help you track, document and report limb volume changes in your patients. The capabilities of LVP6.0 are described in the software using the INFORMATION tab, on the CD itself and also at www.limbvolumes.org.

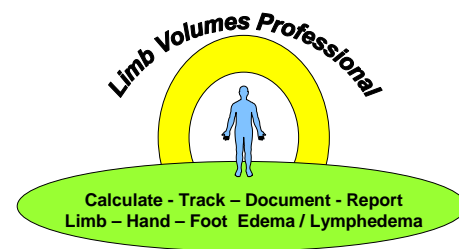

LVP was "field-tested" by experienced therapists and nurses; it got rave reviews. Also, feedback from referring physicians and other health care providers has indicated that the printed/graphic reports of their patient's progress have been extremely useful. They have recommended its use as a "gold-standard" for therapy progress documentation. In fact, LVP has essentially achieved this ranking among therapists worldwide.

INSTALLATION:

To begin installation, put the CD in the drive and wait for it to auto-run. If it does not auto-run then navigate to the CD and double click on the application setupLVP60CD.exe to start the process. During installation LVP6.0 will be put on to your hard drive for future use. Its install location will be C:\Users\Public\LVP60CD. To run LVP6.0 after installation double click on its icon that is placed on your desktop or navigate to the install location and click on LVP60.xls. Note that LVP6.0 uses macros (small bits of code), to run properly. So, if you are asked, please select OK for macros to be enabled. LVP6.0 runs in conjunction with Excel, so you will need to have Excel on your computer. You will be able to run LVP6.0 after you complete the installation. If you are using Excel 2007 then please refer to The information on the CD titled 'UsingLVPwithExcel2007.pdf'.

OPERATION:

Operation is self-explanatory and is displayed in a standard Excel format. Data may be entered, saved by any file name and printed exactly as you would do in Excel. A tutorial and various instructional materials and references are located on the CD and are also available at www.limbvolumes.org. These should be very helpful.

SAVING AND RECALLING

After entering patient data into LVP60.xls you should save the file using a distinct patient identifier of your choosing – such as Mary Smith or MS123456 etc. This is done using the Save As function in the Excel menu. Do not use the save function when working in LVP50.xls as this will write the data you have entered into this master file (LVP60.xls) which you want to keep free of data for ease of subsequent use. However, when recalling a patient's file, with previously entered data, this patient file may be re-saved with the save function.

REGISTRATION and ACTIVATION

This is a single machine license version of LVP6.0, so please install it on the machine you will use. Once installed, you can use it 12 times before it must be permanently activated. During installation a serial number is generated and displayed and is also placed in the HELLO page of the installed software. You will need it and the CD# that appears on the front of the CD for registration. You can register at www.limbvolumes.org. Your permanent activation code will be sent to the e-mail address supplied during registration. When the 12 uses have expired the software automatically asks you to enter an activation code. This is when you should enter it in the box presented. You can activate sooner by using the ACTIVATE tab within the software.

There is a large amount of descriptive material that can be accessed using the INFORMATION tab within the software and at www.limbvolumes.org. If you have any problems, comments or questions please do not hesitate to contact us via e-mail at support@limbvolumes.org. We recognize that any product can be improved to meet specific needs of our clients and we would welcome any feedback that you may wish to share.